

## Certificate of Analysis

### CORNING® COLLAGEN I High Concentration (HC), RAT TAIL

Collagen I is found in most tissues and organs, but is most plentiful in dermis, tendon and bones. Type I molecule is a heterotrimer [ $\alpha_1(I)_2 \alpha_2(I)$ ] of 300 nm length being composed of two  $\alpha_1(I)$  chains and one  $\alpha_2(I)$  chain.<sup>1</sup> Collagen binding integrin receptors are  $\alpha_1 \text{Beta}_1$ ,  $\alpha_2 \text{Beta}_1$ , and  $\alpha_3 \text{Beta}_1$ .<sup>2</sup> When used as a gel, collagen facilitates successful adaptation *in vitro* culture and enhances expression of cell-specific morphology and function. The Collagen HC is used in applications where a sturdy gel provides maximal support to maintain the three dimensional structure.<sup>3-5</sup>

CATALOG NUMBER:	354249	LOT NUMBER:	9238010
SOURCE:	Rat tail tendon		
QUANTITY:	100 milligrams protein (measured by Pyrochemiluminescence)		
CONCENTRATION:	8.63 mg/mL		
FORMULATION:	0.02 N Acetic Acid		
USE:	Corning Collagen I HC, rat tail, is used as a three dimensional gel. Please see reverse for gelling procedures. Use these as guidelines only - we recommend that each investigator empirically determine the optimal conditions for their unique applications.		
QUALITY CONTROL:	<p>≥ 90 % by SDS PAGE.</p> <p>This product has been tested for its ability to promote the attachment and spreading of HT-1080 Human Fibrosarcoma cells.</p> <p>Corning Collagen I HC, rat tail, is a membrane-filtered (0.2 micron) preparation. Tested and found negative for the presence of bacteria, fungi and mycoplasma.</p>		
STORAGE:	<p>Stable when stored at 2-8°C. <b>DO NOT FREEZE.</b></p> <p>On release this product has been successfully gelled over a wide range of dilutions and will form a gel up to a concentration of 0.3 mg/mL. Further dilution may decrease the rigidity of the gel as will the time from manufacture.</p>		
EXPIRATION DATE:	October 26, 2021		
REFERENCES:	<ol style="list-style-type: none"> <li>1. Linsenmayer, T.F., Collagen, in Cell Biology of Extracellular Matrix (ed., E.D. Hay) pp 5-37, Plenum Press, NY (1991).</li> <li>2. Chan, B.M., and Hemler, M.E., <i>J. Cell Biol.</i>, <b>120</b>:537 (1993).</li> <li>3. Gautreau, A., et.al., <i>PNAS</i>, <b>96</b>:7300 (1999).</li> <li>4. Abir, R., et.al., <i>Hum. Reprod.</i>, <b>14</b>:1299 (1999).</li> <li>5. Abir, R., et.al., <i>Fertil. Steril.</i>, <b>75</b>:141 (2001).</li> </ol>		
SAFETY RECOMMENDATION:	Handle in accordance with good industrial hygiene and laboratory safety practices.		

### Suggested Coating Procedures

Gelling Procedure- Corning Collagen I HC, rat tail, will gel when its pH is brought to alkalinity using the procedure below. Please use this as a guideline for determining the optimum concentration for your application.

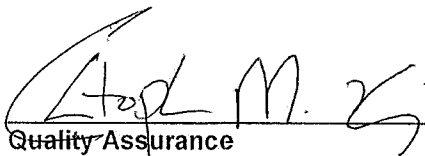
- 1) Prepare ammonia vapor chamber by taping a sterile 2 inch gauze sponge to the inside lid of a 150 mm petri dish. Saturate the gauze with ammonium hydroxide. Place lid on 150 mm dish and set aside.
- 2) Place an even coating of Corning® Collagen I HC, rat tail, on surface to be coated. Thickness may be varied as desired. 50-100 µl of Corning Collagen I HC, rat tail, is sufficient to coat a 22 mm coverslip. For dishes of 100 mm diameter, add approximately 6.0 mL per dish; for 60 mm dishes add approximately 2.3 mL and for 35 mm dishes add approximately 1.0 mL.
- 3) Transfer coated coverslips or dishes with lids off to ammonia vapor chamber and expose for three minutes.
- 4) Soak coated coverslip or dishes in sterile dH<sub>2</sub>O for 30 minutes (5 mL for 35 mm dishes, 10 mL for 60 mm dishes, etc.). Aspirate and replace with 0.5-1.0 mL of sterile dH<sub>2</sub>O and let sit overnight lidded in a laminar flow hood.
- 5) Aspirate the dH<sub>2</sub>O and replace with serum supplemented balanced salt solution and store at 2-8°C.

### Alternate Gelation Procedure for Corning Collagen I HC, Rat tail

- 1.0 Place on ice the following:
  - 1.1 Corning Collagen I HC, rat tail
  - 1.2 Sterile 10X phosphate buffered saline (10X PBS)
  - 1.3 Sterile dH<sub>2</sub>O
  - 1.4 Sterile 1 N NaOH
- 2.0 Determine the final volume of Corning Collagen I HC, rat tail, solution to be used and the desired final collagen concentration.
- 3.0 Place on ice a sterile tube of sufficient capacity to contain the final volume of Corning Collagen I HC, rat tail.
- 4.0 Perform the following steps using aseptic technique in a Class 100 Hood.
  - 4.1 Add to the tube the following volume of 10X PBS:
 
$$\frac{\text{Final Volume}}{10} = \text{mL 10X PBS}$$
  - 4.2 Calculate the volume of Corning Collagen I HC, rat tail, to be used (do not add to the tube until step 4.6):
 
$$\frac{\text{Final volume} \times \text{Final collagen concentration in mg/mL}}{\text{Concentration in bottle (see lot specific spec. sheet)}} = \text{volume collagen to be added}$$
  - 4.3 Add to the 10X PBS the following volume of sterile ice cold 1 N NaOH:
 
$$(\text{volume collagen to be added}) \times 0.023 \text{ mL} = \text{volume 1 N NaOH}$$

- 4.4 Add to the 10X PBS/1 N NaOH the following volume of sterile ice-cold dH<sub>2</sub>O:  
 $(\text{Final volume}) - (\text{Volume collagen}) - (\text{Volume 10X PBS}) - (\text{Volume 1 N NaOH}) = \text{Volume dH}_2\text{O to add}$
- 4.5 Mix the contents of tube and hold in ice.
- 4.6 Add the calculated volume of Corning® Collagen I HC, rat tail, and mix. Leave on ice until ready for use.
- 5.0 The Corning Collagen I HC, rat tail, solution can be used immediately or held on ice for 2-3 hours.
- 6.0 When ready to use, aseptically deliver the solution into the cell culture device and allow to gel at 37°C for 30 minutes.

**NOTE:** For more details on Corning Collagen products and technical resources please visit support page at [www.corning.com/lifesciences](http://www.corning.com/lifesciences)

  
Quality Assurance

October 10, 2019  
Date